


## C30 Glassware for Use in Purging Samples

- Large selection of round bottom flasks and Erlenmeyer flasks
- For use in purging large volume liquid or solid samples
- Adapters available for trapping on S.I.S. GLT desorption tubes
- Constructed of Borosilicate glass
- Flask necks are 24/40 ground glass female fittings



The following glassware can be used in conjunction with the S.I.S. purge head adapter to permit the use of the S.I.S. tube style liquid purge head for purging liquid and solid samples and trapping volatiles purged from the samples onto the S.I.S. GLT Thermal Desorption Tubes. The following glassware all have 24/40 ground glass necks. The PTFE Purge Head Adapter seals into this ground glass fitting and its 1/4" diameter 'O' Ring seal accepts the Tube Style Liquid Purge Head. This Tube Style Purge Head described previously permits both the introduction of the purge gas to purge the volatiles from the sample and the trapping of these volatiles on the GLT desorption tubes.

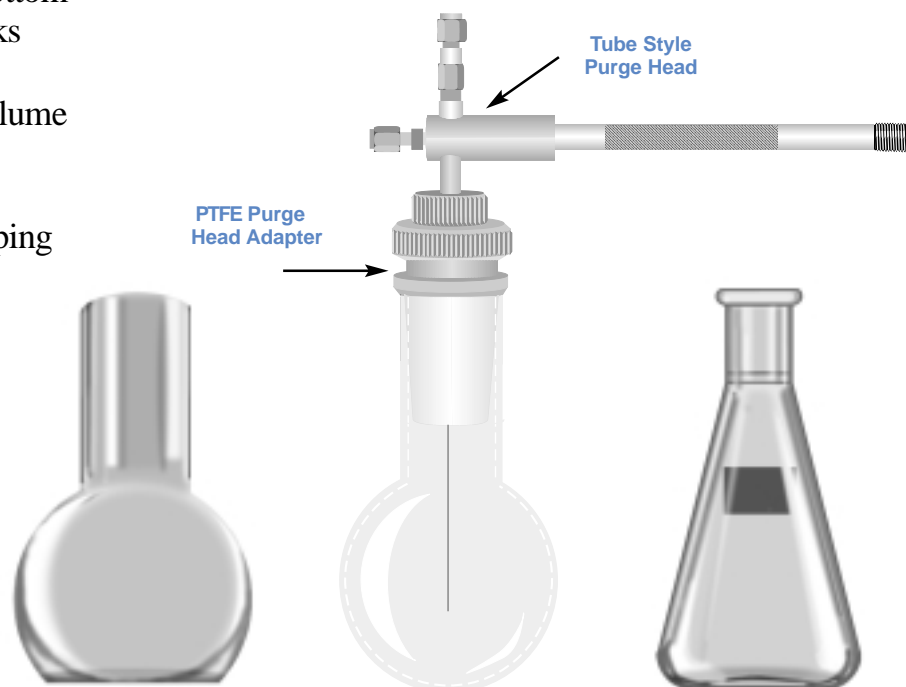

### Round Bottom Flasks

Part No.	Description	Size (ml)	Price
LG7220140	Round Bottom Flask	50 ml	
LG7220142	Round Bottom Flask	100 ml	
LG7220146	Round Bottom Flask	200 ml	
LG7220148	Round Bottom Flask	250 ml	
LG7220152	Round Bottom Flask	500 ml	
LG7220154	Round Bottom Flask	1000 ml	
LG7220156	Round Bottom Flask	2000 ml	

### Flat Bottom Flasks

428126	Flat Bottom Flask	125 ml	
428146	Flat Bottom Flask	250 ml	
428164	Flat Bottom Flask	500 ml	

### Erlenmeyer Flasks, Flat Bottom

635926	Erlenmeyer Flask	125 ml	
LG7750130	Erlenmeyer Flask	500 ml	

### Adapters and Accessories

164372	Purge Head Adapter, PTFE 24/40 Liquid Purge Head, Tube style		
--------	--	--	--