


D22 GC Accessories

SGE Capillary washing/coating reservoir

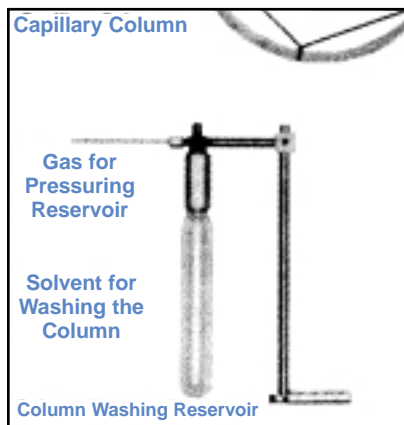

- Ideal for washing and coating capillary columns.
- Removes column contamination through solvent washing
- Used to regenerate a contaminated column

The 10mL Capillary Washing/Coating Reservoir is designed to push solutions through a capillary column. The Reservoir is ideal for washing columns to remove contamination and for custom capillary phase coating.

One end of the column is inserted into the reservoir and sealed using a graphite ferrule GFF/008. Gas pressure, applied through a side arm is used to force the solution through the capillary.

Part No.	Order code	Description	Price
0625026	CWR-10	Washing and Coating Reservoir	
Replacement Parts			
0418490		TCS Pre-drilled Septa (4mm)	
072650		PTFE 1/16" Sealing Ring	

Soap Bubble Flow Meters



- useful for measuring any gas flow rate.
- ideal for measuring GC split and septum purge flow rates

Part No.	Description	Capacity	Price
062505	Bubble Flow Meter	100µl	
062506	Bubble Flow Meter	500µl	
062503	Bubble Flow Meter	50ml	
0625031	Rubber Bulb for Soap bubble Meter	50ml	

Gas & Vapor Field Sampler

- gas and vapor collection
- ideal for field analysis

The Gas and Vapor Field Sampler enables spot collection of environmental gas sampled or low pressure gas streams in 200 ml impervious resealable bags. The undiluted sample can then be simply raised to a required pressure by pressurizing the outside of the bag in the canister.

Part No.	Description	Qty	Price
090111	Gas & Vapor Sampling Kit		
Replacement Parts			
090112	Sample Bag	10	
090113	Septa	5	

