


CDS Pyroprobe Repair

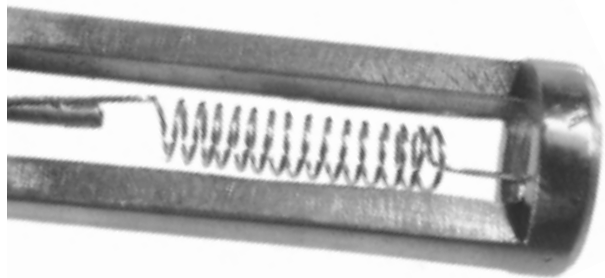

Both the coil and ribbon style pyrolysis probes for the CDS pyroprobe system can be repaired by S.I.S.™ The repair includes the removal of the old heater coil or ribbon, thorough cleaning by sandblasting of the probe tip, polishing the probe rod, replacing the platinum ribbon or coil, and checking for electrical continuity. Repaired probes are not calibrated.

Part No.	Description	Price
PPR-1-C	Pyrolysis Probe Repair with Coil	
PPR-1-R	Pyrolysis Probe Repair with Ribbon	

Replacement Platinum Coils



S.I.S. now supplies platinum replacement coils for the CDS pyroprobe. Coil is made from 0.014" diameter platinum wire.

Part No.	Description	Price ea.
CDS-DC	CDS Pyroprobe Coil	

Platinum Wire and Ribbon for CDS Pyroprobe

For those users who have their own spotwelder, we offer the platinum wire and ribbon from which you can make your own heater coils or ribbons for the CDS pyroprobe.

Part No.	Description	Price/ft
W214	.014" dia. platinum wire	
W215	.0005" x .090" platinum ribbon	

Quartz Tubes for Pyrolysis Probes

There are two styles of quartz tubes. Tubes with both ends open or with one end sealed are available. Both are constructed of high purity tubing. These tubes replace CDS#0100-1060 and have the following dimensions: O.D. - 2.5mm I.D. - 1.9mm Length - 25mm.

The standard quartz tube for the CDS pyrolysis coil probe is open on both ends. In this tube, the sample to be analyzed is held in place with quartz wool. As a second option, we are manufacturing these same quartz tubes with the bottom end of the tube sealed. These tubes can be used in applications where powders or liquids are analyzed where it is not possible to sandwich the sample with quartz wool.

Other sizes and types of quartz and glass vials are available. See the sample vial pages in this catalog for complete details.



Part No.	Description	Price
CDS-QV-20	20 unsealed quartz tubes	
CDS-QV-100	100 unsealed quartz tubes	
CDS-QS-20	20 sealed quartz tubes	
CDS-QS-100	100 sealed quartz tubes	

For Fast and Accurate Delivery
FAX your Orders.....908-806-6631
or order online at:
www.sisweb.com

